

# Educating the Educator: Use of Pulse Oximetry in Athletic Training

David C. Berry, PhD, ATC\*, S. Robert Seitz, MEd, RN, NREMT-P†

\*Saginaw Valley State University, University Center, MI, †University of Pittsburgh, Pittsburgh, PA

## INTRODUCTION

The 5th edition of the Athletic Training Education Competencies expanded the scope of knowledge and skill set of entry-level athletic trainers related to the domain of Acute Care of Injuries and Illnesses.<sup>1</sup> One of these major changes includes the introduction of adjunct airway techniques,<sup>1,2</sup> such as oropharyngeal and nasopharyngeal airways and suction. Another new addition includes the use of pulse oximetry to assess oxygen saturation levels (Table 1). The addition of these acute care knowledge and skill competencies were incorporated into the athletic training curriculum to reflect the National Athletic Trainers' Association's currently published position statement recommendations, as well as those in development.<sup>1</sup> The purpose of this column is to provide athletic training educators (ATE) with a resource on how to use and teach pulse oximetry within their educational programs.

## WHAT IS PULSE OXIMETRY?

Pulse oximetry in the pre-hospital setting provides a rapid, noninvasive measurement and monitoring technique to estimate peripheral oxygen saturation (SpO<sub>2</sub>) of hemoglobin,<sup>3-5</sup> and provides a measure of cardio-respiratory function.<sup>6,7</sup> The use of pulse oximetry also establishes the effectiveness of medical interventions<sup>8-10</sup> through identifiable improvements in baseline SpO<sub>2</sub> levels and can quickly detect hypoxemia<sup>6</sup> (ie, insufficient oxygenation of arterial blood). However, "unlike the thermometer, pulse oximeters are complex machines and the physiological processes they assess are also complex."<sup>4</sup> (p. 24) Thus, it is necessary to understand the underlying physiology of oxygenation that affects arterial oxygen saturation in order to properly interpret the results and guide clinical decision-making. Athletic trainers must understand that oxygen saturation values obtained from pulse oximetry are only ONE component of a complete assessment of the patient's oxygenation status.<sup>5</sup>

## BLOOD GASES AND OXYGENATION

### Oxygen and Carbon Dioxide

Blood gases assist with determining oxygenation, or the exchange of oxygen (O<sub>2</sub>) molecules (moving into the blood from the lungs) and carbon dioxide (CO<sub>2</sub>) molecules (moving into the lungs from the blood), and acid/base (pH) balance. The measurement of blood gases is important in cases of poor oxygenation or pH imbalances where patients demonstrate dyspnea, apnea, nausea or vomiting. These measurements are obtained either noninvasively, with the use of a pulse oximeter to estimate only peripheral oxygen saturation, or invasively (outside the scope of practice of an athletic trainer), with an arterial blood gas measurement, an authentic measurement of respiratory function.

Table 1. Acute Care of Injuries and Illnesses Competencies focusing on Pulse Oximetry

Number	Competency
AC-6	When appropriate, obtain and monitor signs of basic body functions including pulse, blood pressure, respiration, pulse oximetry, pain and core temperature. Relate changes in vital signs to the patient's status.
AC-18	Access oxygen saturation using a pulse oximeter and interpret the results to guide decision-making.

Used with permission from the National Athletic Trainers' Association

Dr. Berry is an associate professor and the professional athletic training education program director at Saginaw Valley State University. Please address all correspondence to David C. Berry, PhD, ATC, Saginaw Valley State University, 7400 Bay Road, University Center, MI 48710. [dcberry@svsu.edu](mailto:dcberry@svsu.edu)

## Full Citation:

Berry DC, Seitz RS. Educating the educator: use of pulse oximetry in athletic training. *Athl Train Educ J*; 2012;7(2):74-80.

Several specific blood gas measurements are possible when drawing arterial blood, making it the gold standard in determining oxygenation. These measures are interrelated with other component tests of blood gases and the results must be considered together when assessing the overall status of a patient. The first value, Partial Pressure of Oxygen In Arterial Blood, or  $\text{PaO}_2$ , is used to assess the pressure of oxygen dissolved in blood and the ability of oxygen to move from the alveoli sac of the lungs into the blood (ie, measurement of gas exchange). Second, Partial Pressure of Carbon Dioxide, or  $\text{PaCO}_2$ , measures the amount of carbon dioxide dissolved in the blood and how well carbon dioxide moves out of the body (ie, assesses adequacy of ventilations). Third, Arterial Oxygen Saturation, or  $\text{SaO}_2$  (oxygen saturation), is a direct measure of the percent of oxygen bound to hemoglobin via a blood draw. Finally,  $\text{SpO}_2$  is the estimated oxygen saturation of hemoglobin as measured by a pulse oximeter and given as a percentage, with normal around 96%. The relationship of  $\text{PaO}_2$ ,  $\text{PaCO}_2$  and  $\text{SaO}_2$  is outlined in Table 2. Several other blood gases values do exist, but are beyond the scope of this article.

### Oxygen Dissociation and Its Effects on Pulse Oximetry

As a foundation for utilizing pulse oximetry, a review of tissue perfusion physiology is also warranted. In addition to an adequate availability of oxygen, 3 other factors impact pulmonary perfusion: (1) an adequate blood volume, (2) intact pulmonary capillaries, and (3) the efficient pumping of blood by the heart.<sup>11</sup> Oxygen is transported in the blood by red blood cells bound to hemoglobin (oxyhemoglobin or  $\text{HbO}_2$ ), or is dissolved in the plasma (deoxyhemoglobin), 98% and 2% respectively.

Oxygen dissociation describes the equilibrium of oxyhemoglobin and deoxyhemoglobin. In areas containing a high partial pressure of  $\text{O}_2$  (ie, lungs),  $\text{O}_2$  binds to hemoglobin to form oxyhemoglobin. The partial pressure of  $\text{O}_2$  will decrease as oxyhemoglobin moves into body tissues that are deprived of oxygen, releasing  $\text{O}_2$  for metabolic use and acquiring  $\text{CO}_2$  for transport back to the alveoli sac. While a strong positive correlation between the percentage of available hemoglobin that is bound to oxygen ( $\text{SaO}_2$ ) and  $\text{SpO}_2$  exists,<sup>12</sup> the assurance of adequate ventilation (carbon dioxide and acid base balance) must be continually monitored in the acutely ill or injured patient in order to avoid overestimation of  $\text{SpO}_2$ , especially in the presence of hypoxemia.<sup>13</sup>

In the presence of increased blood carbon dioxide ( $\text{PaCO}_2$ ), hemoglobin will lose its affinity for oxygen. As  $\text{CO}_2$  levels increase in the tissues, more  $\text{O}_2$  is released from hemoglobin, thus decreasing oxyhemoglobin levels and lowering the pulse oximeter reading (right shift on the oxygen dissociation curve). A loss of affinity for oxygen is also noted with an increase in body temperature, an increase in 2,3-diphosphoglycerate (2,3-DPG), an organic phosphate salt found in red blood cells, or a decrease of blood pH (hydrogen ion concentration). A shift in the oxygen dissociation curve to the left indicates opposite conditions exist. Figure 1 represents the relationship between the partial pressure of  $\text{O}_2$  and oxyhemoglobin, and Table 3 outlines the interpretation for  $\text{SpO}_2$  readings.

Table 2. Normal Blood Gas Values

Gas	Definition	Normal Values
$\text{PaO}_2$	Partial pressure of oxygen that is dissolved in arterial blood	80 to 100 mm Hg
$\text{PaCO}_2$	Amount of carbon dioxide dissolved in arterial blood	35-45 mm Hg
$\text{SaO}_2$	Arterial oxygen saturation	95% to 100%
$\text{SpO}_2$	Estimated of arterial oxygen saturation	95% to 100%

These are average reference ranges. Values may vary between analyzers, laboratories, and reference manuals.

Table 3. Interpretation of Pulse Oximetry Readings.<sup>11</sup>

$\text{SpO}_2$ Reading %	Interpretation
95-100	Normal
91-94	Mild Hypoxemia
86-90	Moderate Hypoxemia
<85	Severe Hypoxemia

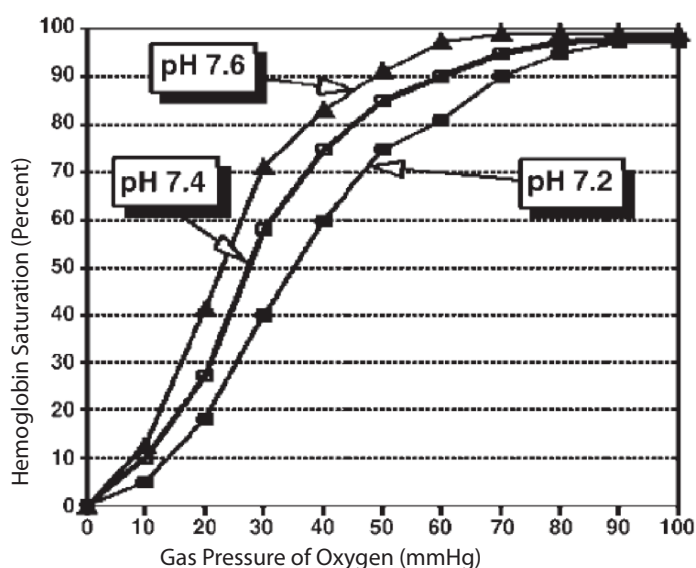


Figure 1. Oxygen Dissociation Curve.<sup>11</sup>

Note: A shift in the curve to the right indicates a decrease in available oxyhemoglobin. Reproduced, with kind permission, from Dr. David McShaffrey, Marietta College, Marietta, OH <http://www.marietta.edu/~mcshaffd/index.html>

Unit Design

Developed in the mid-1970s, pulse oximeters are small hand-held devices that indirectly measure (estimate) oxygen saturation levels and heart rate. The pulse oximeter functions based on the principle that saturated hemoglobin and desaturated hemoglobin absorb light at different frequencies.<sup>4,5,14</sup> Pulse oximeter units use light-emitting diodes containing two light sources, visible (red) and infrared light, and a photo sensor. The saturated hemoglobin absorbs visible light, while desaturated hemoglobin absorbs the infrared light. The light source and the photo sensor are placed on either side of an artery (eg, radialis indicis artery of the index finger). When light from the light-emitting diodes is transmitted through an artery's pulsating blood flow (normally for 5 pulses) and surrounding tissues, the hemoglobin (saturated and desaturated) absorbs the different light frequencies.<sup>14</sup> The difference in the amount of light absorbed by the saturated and desaturated hemoglobin, as detected by the photo sensor, is then converted into a digital value representing the percentage of hemoglobin saturated with oxygen. This conversion is based on a complex algorithm for calculating exact saturation and separating arterial pulsation from motion artifact,<sup>14</sup> and is a close estimation to the arterial oxygen saturation in the hemoglobin.<sup>15</sup>

Pulse Oximetry Measurement Accuracy

Non-physiological Factors

The accuracy of SpO<sub>2</sub> measurements can be affected by both non-physiological and physiological variables. Two non-physiological variables an athletic trainer (AT) can control are frequency of unit calibration and ensuring proper placement of the unit on the finger or earlobe. New unit measurement accuracy is approximately ±2% for values over 90% and higher when values are reported below 80% saturation.<sup>4</sup> To maintain a high level of accuracy, units should be calibrated on a yearly basis. When a unit is placed improperly on the finger (or other body part being assessed) the light does not pass through the vascular bed correctly when traveling from the light-emitting diode to the photo sensor.<sup>15</sup>

Table 4. Factors Affecting the Placement and Function of a Finger Tip Pulse Oximeter

Trauma to the digits (ie, contusion)
Decreased peripheral pulses
Decreased core temperature
Decreased blood pressure
Exposure to excessive ambient light (eg, too large of a sensor over a finger)
Excessive movement
Presence of nail polish
Incorrect sensor size (eg, too large sensor will not stay in place, too small of a sensor may constrict arterial flow and provide false low readings).

Because pulse oximeters use light and color to determine oxygen saturation levels<sup>4</sup> the amount of bright ambient light seen by the photo sensor, a dirty sensor,<sup>7</sup> and dried blood can affect the accuracy of the SpO<sub>2</sub> measurement. Patient movement or dislodging the unit (known as motion artifact) experienced with a seizure, hypothermia (ie, shivering), and cardiopulmonary resuscitation also reduces the accuracy of the unit.<sup>15,16</sup> Therefore, when using a pulse oximeter, particularly a figure tip unit, it is necessary to consider when and where to place the unit for the most accurate reading, (Table 4) as a poor signal will result in inaccurate or absent readings.

The effects of nail polish<sup>17-23</sup> and acrylic nails<sup>24</sup> on pulse oximetry have also been studied extensively. Early studies (1980s) on non-hypoxic patients found clinically significant decreases in SpO<sub>2</sub> values<sup>19,21,22</sup> particularly with the nail polish colors of blue, black, green, and purple. Later studies (2000s) examining the effects of nail polish color demonstrate more accurate SpO<sub>2</sub> measurements, with SpO<sub>2</sub> changes less than 2%, compared to the control groups when black, blue, and purple nail-polish were used.<sup>17,18,20</sup> In the mildly hypoxic patient (induced by changes in altitude) nail-polish color (black, blue, purples, green, brown/maroon, white glitter, red, orange, and pink) appear to have no significant effect on SpO<sub>2</sub> values.<sup>23</sup> In the presence of artificial nails there appears to be no significant difference in SpO<sub>2</sub> values between the unpolished acrylic fingernails (97.53%) and natural unpolished fingernails (97.33%).<sup>24</sup>

Physiological Factors

One significant physiological factor affecting the pulse oximeter's measure of cardiorespiratory function is a low perfusion state. Because SpO<sub>2</sub> measurements are a function of the artery's pulsating blood flow, a patient must have sufficient perfusion in the monitored area (eg, index finger). Patients with a weak or absent peripheral pulse, such as those with hypotension, hypovolemia, hypothermia, or cardiac arrhythmia may present with inaccurate measurements<sup>15</sup> due to a lack of an adequate peripheral pulse. Additionally, when using an automatic blood pressure cuff a pulse oximeter should not be placed ipsilaterally because blood flow to the finger will be greatly reduced whenever the cuff inflates.<sup>15,25</sup>

Patients who smoke and/or who are diagnosed with anemia may have falsely high SpO<sub>2</sub> measurements.<sup>4,15</sup> Carbon monoxide, a by-product of smoking, is a colorless and odorless gas that competes for the same binding sites with oxygen on hemoglobin. In fact, the carbon monoxide bind to hemoglobin is approximately 250 times stronger than that of an oxygen molecule.<sup>15</sup> When carbon monoxide binds to hemoglobin it alters the color of hemoglobin making the pulse oximeter unable to accurately detect hemoglobin carrying oxygen and carbon monoxide.<sup>4,15</sup> Anemia affects pulse oximeter measurements in patients during hypoxemia. The causation is not completely understood, but is believed to "be due to photon scattering of light and a shift in red light wavelength increasing its absorption."<sup>26</sup> (p. 43) Other physiological factors affecting SpO<sub>2</sub> measurements can be found in Table 5.

Indications, Contraindications, and Precautions

The use of pulse oximetry does not preclude the necessity of a thorough comprehensive and ongoing assessment of your patient. A clinical assessment of adequate oxygenation and perfusion status is always required, with treatment based on clinical findings and patient history. Waiting for the pulse oximeter to indicate less

Table 5. Physiological Factors Affecting SpO<sub>2</sub> Measurements

Arterial blood flow to the vascular bed
Hemoglobin level
Temperature of the digit or the area where the oximetry sensor is located
Patients's oxygenation ability
Percentage of inspired oxygen
Venous return at the probe location
Presence of anemia, carbon monoxide, cyanide or increased CO <sub>2</sub> (carbaminohemoglobin)
Abnormal pH
Increase or decrease in 2, 3-DPH, CO <sub>2</sub> or core body temperature

than optimal tissue oxygenation places you, and the patient, in a reactive versus proactive response situation. As indicated previously, the assurance of adequate ventilation remains a priority regardless of the value obtained by pulse oximetry.

Indications for the use of pulse oximetry include the need to: (1) monitor the adequacy of arterial oxyhemoglobin saturation in any individual who presents with a complaint or assessment findings related to altered mental status, respiratory distress or signs and symptoms of shock, (2) quantify the response of arterial oxyhemoglobin saturation to therapeutic intervention, (3) quantify the response of arterial oxyhemoglobin saturation during a diagnostic procedure, and (4) comply with institutional or facility policies and procedures.<sup>27</sup> Other than the potential for a formation of pressure ulcers secondary to long-term use, no relative contraindications exist.<sup>27</sup>

Even though pulse oximetry is a safe procedure, false-negative results for hypoxemia and/or false-positive results for normoxemia or hyperoxemia can lead to inappropriate patient treatments.<sup>27</sup> Conditions where hemoglobin is bound by carbon monoxide, cyanide, or CO<sub>2</sub> (carbaminohemoglobin), may present with tissue hypoxia even with a high measurement reading by the pulse oximeter. Ensure worn probes are replaced with new probes designed specifically for your unit and that the unit is calibrated at least yearly. Injury from an electrical shock and/or burns can occur when incompatible probes are used.<sup>27</sup>

Application

The steps describing how to properly use a fingertip pulse oximeter are found in Table 6. Once placed on the finger the AT should monitor the position and movement of the pulse oximeter. Excessive movement/shaking of unit can result in unreliable measurements as well as exposure to excessive bright light, which can result in inflated values. If either of these events occurs, consider moving the unit to a different site and shielding the unit with a blanket or towel, respectively. The following information should be documented during the course of care: (1) patient's SpO<sub>2</sub> level and location the reading was obtained (finger, earlobe), (2) patient's heart rate, (3) patient's assessment at the time of measurement, and (4) any other factor affecting appropriate assessment. At the conclusion of

care, remove the unit and clean according to the manufacturer's directions.

EDUCATIONAL CONSIDERATIONS

Equipment

Today there are many different makes and models of pulse oximeters available in the marketplace. Most models provide a visual (bar graph, waveform) digital output and audible display of arterial pulsations and heart rate<sup>15</sup> particularly when the SpO<sub>2</sub> or heart rate falls below a predefined level. The size and weight of the pulse oximeter should be considered, as lighter units normally result in less unit movement of the unit and therefore more accurate results. Some models offers resistance to light interference to increase the sensitivity of the oxygen saturation measurement and many manufacturers equip pulse oximeters with a signal quality monitor and alarm.<sup>26</sup> Some models offer software to upload and analyze up to 24 hours of data at a time, while others offer the ability to transmit measurement with Bluetooth wireless technology. The selection of a particular unit will depend on the intended use of the unit and setting in which it is used.<sup>15</sup>

Practicing pulse oximetry in the pre-hospital or educational setting requires very little equipment – normally just a fingertip pulse oximeter and a carrying case to protect the unit. Units can range in price from \$69.00 to \$385.00. Price is obviously contingent on the features attached to the unit as well as the quality of the light-emitting diodes and photo sensors (Table 7). An Internet search for “finger pulse oximeter” using Google Shopping results in units ranging in price from \$22.00-\$382.00. The moral of the story here is, “Shoppers beware!”

Teaching the Skill

A tremendous benefit to utilizing this technology resource is the ease of education. From a cognitive aspect, just a couple of questions on what is pulse oximetry, how it works, and how to interpret the information may be easily developed. The psychomotor education component is accomplished through actual application of the device in the classroom and clinical setting. The non-invasive nature along with the ability to alter information output through restriction of blood flow (eg, applied to an extremity while taking a blood pressure) easily demonstrates varying feedback of audio and waveform output as well as digital readouts. Remember, there may be restrictions to utilizing educational equipment in the actual clinical setting dependent upon organizational structure and student policy and procedures. Please refer to your individual biomedical department or administrative resource to determine specific restrictions, if any.

CONCLUSION

The concepts of pulse oximetry, especially those required by educational competencies, should not be viewed as a challenge; rather they should be viewed as an advancement of the profession. Yes, it does add to already long list of knowledge and skills required for an entry-level AT, however, these are also skills required for proper emergency care of a patient. If ATs are truly going to be the “first” on the scene of an emergency situation, then having the requisite training on one of the most fundamental skills taught in every emergency medical responder course is important.

Table 6. Using a Fingertip Pulse Oximeter



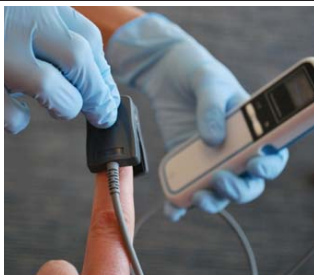

Steps	Procedure	Image
1.	Assume body substance isolation (BSI).	
2.	Select the appropriate size fingertip pulse oximeter and sensor to optimize signal capture and minimizes artifact-related difficulties.	
3.	Turn on the pulse oximeter and ensure proper functioning of the unit.	
4.	Select the desired sensor site; normally the index finger of hand without a sphygmomanometer attached.	
5.	Apply the pulse oximeter in a method that allows the light-emitting diodes to be opposite the photo sensor and shielded from excessive ambient light.  Note: Consider removing fingernail polish as needed. Do not delay immediate care.	
6.	Allow the unit to measure the pulsating blood flow (normally 5 pulses).	
7.	Determine accuracy of by comparing the manually measured numeric heart rate value to the unit's monitored heart rate.	
8.	Document the baseline results.	
9.	Monitor changes in the patient pulse oximetry values and respiratory status.  Do not wait for less-than-optimal tissue oxygenation as reported via the pulse oximeter before providing medical interventions.	

Table 7. Pulse Oximeter Finger Unit Resources.

Company	Web Address	Equipment	Retail Price
Savelives.com	<a href="http://www.savelives.com">http://www.savelives.com</a>	SPO Medical PulseOx 5500 Finger Unit	\$189.00
		SPO5500 Carry Case	\$15.00
		BCI Digit Pulse Oximeter	\$380.00
		BCI Fingerprint Pulse Oximeter	\$665.00
		BCI Fingerprint Nylon Carry Case	\$35.00
		Model 9500 Finger Sized Pulse Oximeter	\$385.00
		Model 9500 Belt Loop Carry Case	\$17.50
eNasco	<a href="http://www.enasco.com">http://www.enasco.com</a>	Fingertip Pulse Oximeter	\$88.00
CPR Saver	<a href="http://www.cpr-savers.com">http://www.cpr-savers.com</a>	Finger Pulse Oximeter – Item - MDC22	\$69.95
		Digit® Finger Oximeter	\$319.00
		View 02 Fingertip Pulse Oximeter – Item - 180-5055	\$128.25
		Fingertip Pulse Oximeter - Item - MS-74001	\$225.00
Laerdal	<a href="http://www.laerdal.com">http://www.laerdal.com</a>	Contact Laerdal directly at (877) LAERDAL	

Note: Speak with a company representative to inquire about educational pricing

## REFERENCES

- National Athletic Trainers' Association. Athletic Training Education Competencies. 5th ed. Dallas, TX: National Athletic Trainers' Association; 2011.
- Berry DC, Seitz RS. Educating the educator: teaching airway adjunct techniques in athletic training. *Athl Train Educ J*. 2011;6(2):107-116.
- Brown LH, Manring EA, Kornegay HB, NH P. Can prehospital personnel detect hypoxemia without the aid of pulse oximeters? *Am J Emerg Med*. 1996;14(1):43-46.
- Casey G. Pulse oximetry--what are we really measuring? *Nurs N Z*. Apr 2011;17(3):24-29.
- Schutz SL. Oxygen saturation monitoring by pulse oximetry. In: Lynn-McHale Wiegand DJ, Carlson KK, eds. *AACN Procedure Manual for Critical Care*. Philadelphia, PA: W.B. Saunders; 2001:77-82.
- Saneesh PJ, Kaimar P. Pulse oximeter: a boon or a bane. *Internet J Anesthesiol*. 2009;20(1):15. <http://www.ispub.com/journal/the-internet-journal-of-anesthesiology/volume-20-number-1/article/pulse-oximeter-a-boon-or-a-bane.html>.
- Sims J. Making sense of pulse oximetry and oxygen dissociation curve. *Nurs Times*. 1996;92(1):34-35.
- Aughey K, Hess D, Eitel D, et al. An evaluation of pulse oximetry in prehospital care. *Ann Emerg Med*. 1991;20(8):887-891.
- Grace RF. Pulse oximetry. Gold standard or false sense of security? *Med J Aust*. 1994;160(10):638-644.
- McGuire TJ, Pointer JE. Evaluation of a pulse oximeter in the prehospital setting. *Ann Emerg Med*. Oct 1988;17(10):1058-1062.
- Bledsoe BE, Porter RS, Cherry RA, Snyder SR. Pulmonology. *Paramedic Care: Principles and Practice*. Vol 3. 3rd ed. Upper Saddle River, NJ: Pearson Education, Inc; 2009:15-16.
- Tittle M, Flynn MB. Correlation of pulse oximetry and co-oximetry. *Dimens Crit Care Nurs*. 1997;16(2):88-95.
- Wilson BJ, Cowan HJ, Lord JA, Zuege DJ, Zygun DA. The accuracy of pulse oximetry in emergency department patients with severe sepsis and septic shock: a retrospective cohort study. *BMC Emerg Med*. 2010;10:9.
- Duran DJ, Phillips B, J B. Blood gases: technical aspects and interpretation. In: Goldsmith JP, Karotkins EH, eds. *Assisted Ventilation of the Neonate*. 4th ed. Philadelphia, PA: W.B. Saunders; 2003:297-292.
- Howell M. Pulse oximetry: an audit of nursing and medical staff understanding. *Br J Nurs*. 2002;11(3):191-197.
- Woodrow P. Pulse oximetry. *Nurs Stand*. 1999;13(42):42-46; quiz 47.
- Brand TM, Brand ME, Jay GD. Enamel nail polish does not interfere with pulse oximetry among normoxic volunteers. *J Clin Monit Comput*. 2002;17(2):93-96.
- Chan MM, Chan ED. What is the effect of fingernail polish on pulse oximetry? *Chest*. 2003;123(6):2163-2164.

19. Cote CJ, Goldstein EA, Fuchsman WH, Hoaglin DC. The effect of nail polish on pulse oximetry. *Anesth Analg*. 1988;67(7):683-686.
20. Hinkelbein J, Genzwuerker HV, Sogl R, Fiedler F. Effect of nail polish on oxygen saturation determined by pulse oximetry in critically ill patients. *Resuscitation*. 2007;72(1):82-91.
21. Rubin AS. Nail polish color can affect pulse oximeter saturation. *Anesthesiology*. 1988;68(5):825.
22. White L. Use of pulse oximeter with nail polish. *AORN J*. 1989;50(3):520.
23. Yamamoto LG, Yamamoto JA, Yamamoto JB, Yamamoto BE, Yamamoto PP. Nail polish does not significantly affect pulse oximetry measurements in mildly hypoxic subjects. *Respir Care*. 2008;53(11):1470-1474.
24. Peters SM. The effect of acrylic nails on the measurement of oxygen saturation as determined by pulse oximetry. *AANA J*. 1997;65(4):361-363.
25. Carroll P. Pulse oximetry--at your fingertips. *RN*. 1997;60(2):22-26; quiz 27.
26. Mengelkoch LJ, Martin D, Lawler J. A review of the principles of pulse oximetry and accuracy of pulse oximeter estimates during exercise. *Phys Ther*. 1994;74(1):40-49.
27. American Association for Respiratory Care. AARC (American Association for Respiratory Care) clinical practice guideline. Pulse oximetry. *Respir Care*. 1991;36(12):1406-1409.